

Whether it's viscous, highly mobile or foaming. The pipette for difficult media.



Transferpettor

FIRST CLASS · BRAND

H_2SO_4

Silicone oil

ISOPROPANOL

Blood

Glycerine

Ethanol

Mineral oil

The Transferpettor is ideal for liquids when air displacement pipettes just won't work.

By virtue of its positive displacement principle, the Transferpettor is the right tool for your most-demanding pipetting operations. Convenient and tidy pipetting of difficult media

– with highest precision.

- Viscous
- Highly mobile
- Foaming
- High density
- High vapor pressure



Алматы (7273)495-231
Ангарск (3955)60-70-56
Архангельск (8182)63-90-72
Астрахань (8512)99-46-04
Барнаул (3852)73-04-60
Белгород (4722)40-23-64
Благовещенск (4162)22-76-07
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Владикавказ (8672)28-90-48
Владимир (4922)49-43-18
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89

Россия +7(495)268-04-70

Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Иркутск (395)279-98-46
Казань (843)206-01-48
Калининград (4012)72-03-81
Калуга (4842)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Коломна (4966)23-41-49
Кострома (4942)77-07-48
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Курган (3522)50-90-47
Липецк (4742)52-20-81

Казахстан +7(7172)727-132

Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41
Нижний Новгород (831)429-08-12
Новокузнецк (3843)20-46-81
Ноябрьск (3496)41-32-12
Новосибирск (383)227-86-73
Омск (3812)21-46-40
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Петрозаводск (8142)55-98-37
Псков (8112)59-10-37
Пермь (342)205-81-47


Киргизия +996(312)96-26-47

Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Саранск (8342)22-96-24
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78
Севастополь (8692)22-31-93
Симферополь (3652)67-13-56
Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Сургут (3462)77-98-35
Сыктывкар (8212)25-95-17
Тамбов (4752)50-40-97
Тверь (4822)63-31-35

Тольятти (8482)63-91-07
Томск (3822)98-41-53
Тула (4872)33-79-87
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Улан-Удэ (3012)59-97-51
Уфа (347)229-48-12
Хабаровск (4212)92-98-04
Чебоксары (8352)28-53-07
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Чита (3022)38-34-83
Якутск (4112)23-90-97
Ярославль (4852)69-52-93

This is how it works

Convenient pipetting of difficult media.

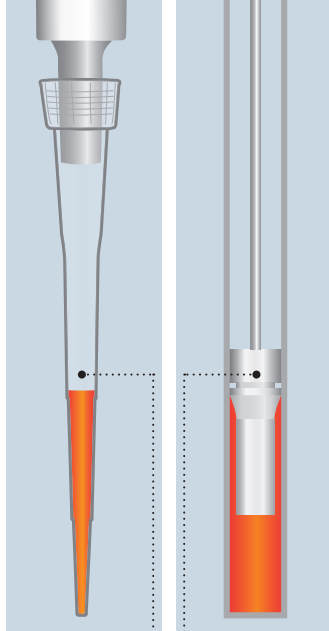

Air-interface principle

Positive displacement principle



Media which tend to foam, such as tenside solutions, pose no problem to the Transferpettor.



Even media with high vapor pressure (up to 500 mbar) such as alcohols, ether and hydrocarbons can be pipetted simply, reliably and with the highest accuracy.



The Transferpettor is ideal for pipetting highly viscous media, such as highly concentrated protein solutions, oils, resins and fats up to a viscosity of 140,000 mm²/s. Media with high density up to 13.6 g/cm³ such as glycerin, mercury, sulfuric acid, etc. can be pipetted without the need to recalibrate.

Transferpettor, Digital-adjustable, DE-M marking

Volume µl	A* ≤ ± % µl	CV* ≤ % µl	Subdiv. µl	Color code	Cat. No.
2.5 - 10	1.0 0.1	0.8 0.08	0.01	orange	7018 07
5 - 25	0.8 0.2	0.5 0.125	0.1	2 x white	7018 12
10 - 50	0.6 0.3	0.4 0.2	0.1	green	7018 17
20 - 100	0.6 0.6	0.4 0.4	0.1	blue	7018 22
100 - 500	0.5 2.5	0.2 1.0	1.0	green	7028 04
200 - 1000	0.5 5.0	0.2 2.0	1.0	yellow	7028 06
1000 - 5000	0.5 25.0	0.2 10.0	10.0	red	7028 10
2000 - 10000	0.5 50.0	0.2 20.0	10.0	orange	7028 12

Transferpettor-Station

The Transferpettor-Station keeps the instruments and accessories ready to hand at the workplace.

The color code clearly identifies which belongs to what, to avoid any confusion.

Transferpettor, Fixed-volume, DE-M marking

Volume µl	A* ≤ ± % µl	CV* ≤ % µl	Color code	Cat. No.
1	4.0 0.04	4.0 0.04	white	7018 42
2	2.5 0.05	2.0 0.04	white	7018 44
5	1.0 0.05	0.8 0.04	white	7018 53
10	1.0 0.1	0.8 0.08	orange	7018 58
20	0.8 0.16	0.5 0.1	black	7018 63
25	0.8 0.2	0.4 0.1	2 x white	7018 64
50	0.6 0.3	0.4 0.2	green	7018 68
100	0.6 0.6	0.4 0.4	blue	7018 73
200	0.5 1.0	0.2 0.4	red	7018 78

Accommodates 4 instruments up to 200 µl with accessories. Pack of 1.

Cat. No. 7019 60



Accommodates 2 instruments 0.5 to 10 ml with accessories. Pack of 1.

Cat. No. 7028 90

* Calibrated to deliver (TD, Ex). Error limits according to the nominal capacity (= maximum volume) indicated on the instrument, obtained with instrument and distilled water at equilibrium with ambient temperature at 20 °C, and with smooth, steady operation. The error limits are within the limits of DIN EN ISO 8655-2. DE-M marking.

A = Accuracy, CV = Coefficient of variation

Алматы (7273)495-231
Ангарск (3955)60-70-56
Архангельск (8182)63-90-72
Астрахань (8512)99-46-04
Барнаул (3852)73-04-60
Белгород (4722)40-23-64
Благовещенск (4162)22-76-07
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Владикавказ (8672)28-90-48
Владимир (4922)49-43-18
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89

Россия +7(495)268-04-70

Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Иркутск (395)279-98-46
Казань (843)206-01-48
Калининград (4012)72-03-81
Калуга (4642)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Коломна (4966)23-41-49
Кострома (4942)77-07-48
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Курган (3522)50-90-47
Липецк (4742)52-20-81

Казахстан +7(7172)727-132

Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41
Нижний Новгород (831)429-08-12
Новокузнецк (3843)20-46-81
Ноябрьск (3496)41-32-12
Новосибирск (383)227-86-73
Омск (3812)21-46-40
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Петрозаводск (8142)55-98-37
Псков (8112)59-10-37
Пермь (342)205-81-47

Киргизия +996(312)96-26-47

Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Саранск (8342)22-96-24
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78
Севастополь (8692)22-31-93
Симферополь (3652)67-13-56
Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Сургут (3462)77-98-35
Сыктывкар (8212)25-95-17
Тамбов (4752)50-40-97
Тверь (4822)63-31-35

Тольятти (8482)63-91-07
Томск (3822)98-41-53
Тула (4872)33-79-87
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Улан-Удэ (3012)59-97-51
Уфа (347)229-48-12
Хабаровск (4212)92-98-04
Чебоксары (8352)28-53-07
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Чита (3022)38-34-83
Якутск (4112)23-90-97
Ярославль (4852)69-52-93